



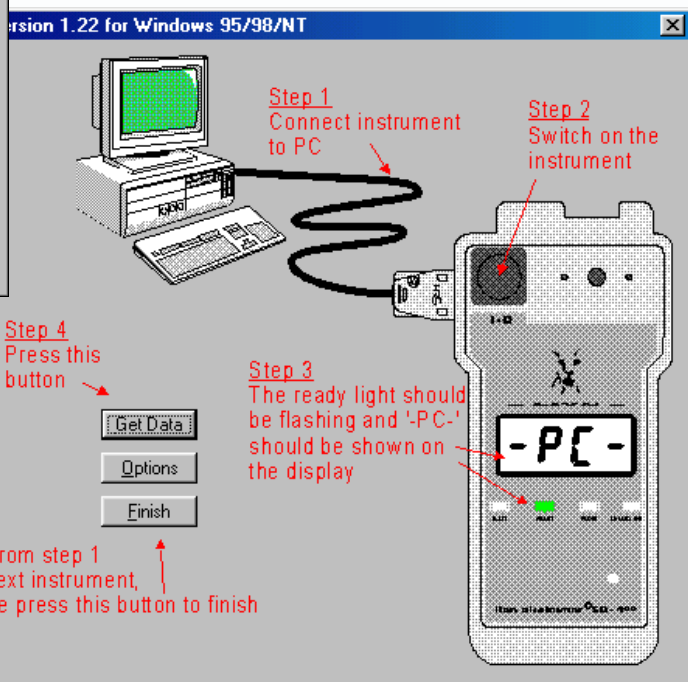

Data-400™

software for the lion alcolmeter® 400

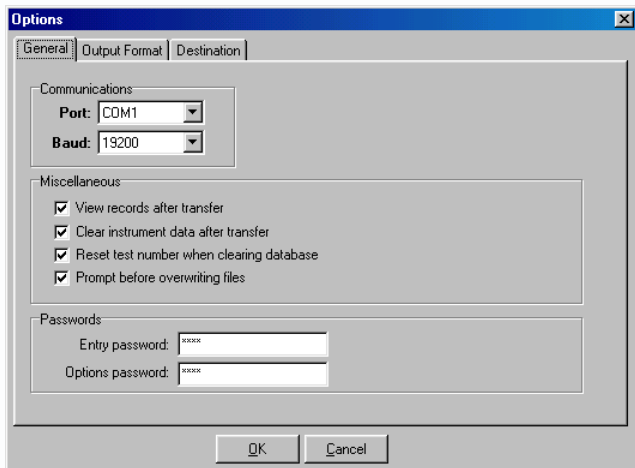

Each time a **lion alcolmeter® 400** tests a subject's breath, **data** is stored in its memory. This data may simply be the alcohol reading, with date and time; or it may, if the instrument has data logging capabilities, include the test location and reason codes, the identity number of the operator, and the sex of the subject as well.

Data-400™ extracts this data from each instrument and stores it on your PC's hard drive. From there it may be printed out or exported to another software package, such as Microsoft Excel* or Lotus 1-2-3*. Data from one instrument may then be merged with data from all others, to build a complete analysis of breath testing in your area.

Technical Information

Windows system requirements: IBM compatible personal computer [minimum 80386X], mouse, at least 3Mb of available hard drive space and 4Mb of available RAM, Microsoft Windows* version 3.1 or higher, VGA or higher resolution display, CD-ROM.

Data output formats: data may be exported using comma, space or Tab separators, according to the requirements of the destination software system.

Supplied with each pack: installation disk, cable for connecting the **400** to the computer, instruction manual.

Important ordering information: when ordering, please specify details of the operating system in your computer, and the software version number in your **400** [displayed on entry into Calibration Check mode].